



Multi-function Holders

68202 Large Probe Holder



Lateral fixation, holds micro syringe needles, large diameter electrodes and flexible drill handles etc, V-groove design with 6mm-12mm diameters.

68106 Headstage Holder



Holder is designed for holding micro amplifier or headstage (Such as TBSI headstage and Plexon microelectrode amplifier). C-shaped head and suitable for diameter from 0.1-7.0mm.

68218 Syringe Holder




Holder is applicable for micro-injection, which holds syringe and needle simultaneously, and the clamping of the syringe is more stable, the syringe holding range is 3mm to 12mm and needle's is 0.3mm to 1.5mm.

68605 Microdrill Holder



Holder secures the microdrill to stereotaxic instruments. Controls the depth of drilling by the manipulator. Easy and accurate operation protects brain tissue from damage caused by excessive drilling.

68397 Level Calibrator - Rat

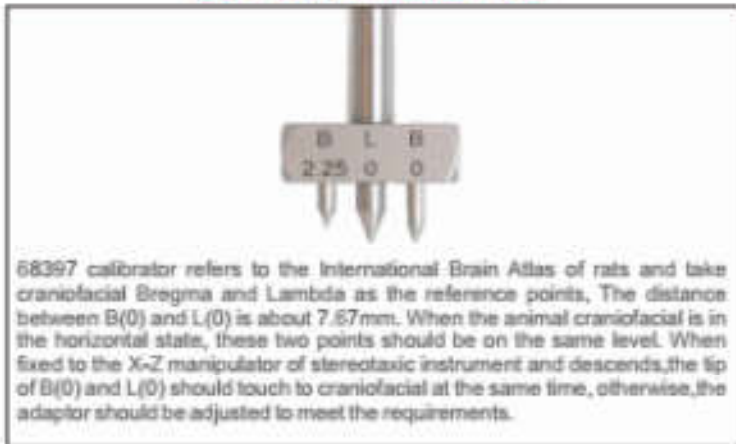

68397 calibrator refers to the International Brain Atlas of rats and take craniofacial Bregma and Lambda as the reference points. The distance between B(0) and L(0) is about 7.67mm. When the animal craniofacial is in the horizontal state, these two points should be on the same level. When fixed to the X-Z manipulator of stereotaxic instrument and descends, the tip of B(0) and L(0) should touch to craniofacial at the same time, otherwise, the adaptor should be adjusted to meet the requirements.